

Support Every Step Of The Way

We're here to guide you every step of the way, offering the support and resources you need to feel confident and informed throughout your journey.



Get In Touch!

Telephone: 01635 887640

Email: service@severnhealthcare.com

Website: www.severnhealthcare.com



Severn Healthcare Technologies Ltd
42 Kingfisher Court,
Hambridge Road,
Newbury, Berkshire
RG14 5SJ

Stay connected with us:



INHEALTH[®]
TECHNOLOGIES



Day&Night[™] HME

INHEALTH[®]
TECHNOLOGIES

Let's Get Started!




Day&Night® HME

It's our newest
innovation that
you can wear
comfortably all
day and all night –
including bedtime,

Love the simplicity.
Love the innovation!


The Day&Night HME cartridge contains a foam filter for heat and moisture exchange, side airflow vents to maximise airflow, and a locking mechanism to secure the device in the open position to prevent occlusion as desired.

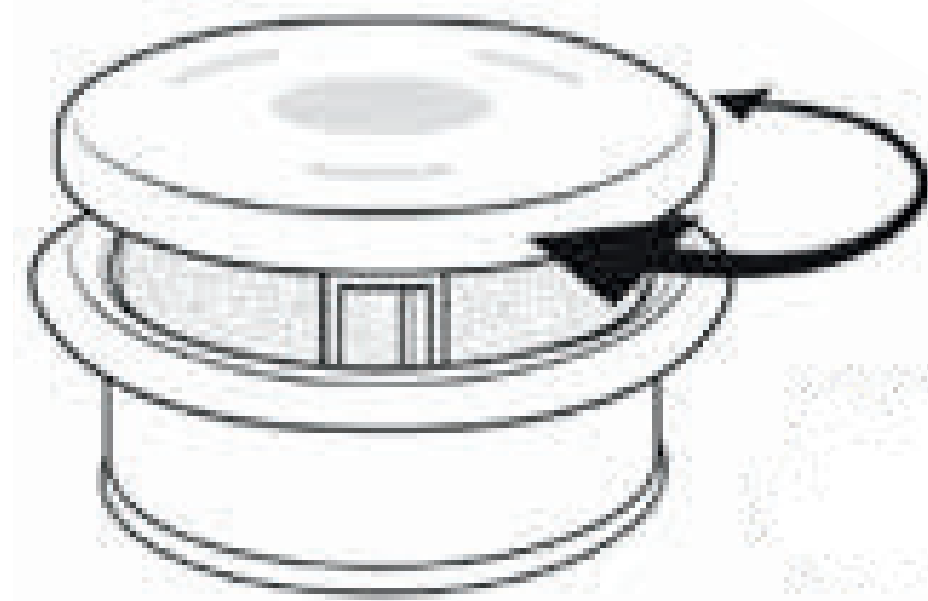

STEP 1:

Connect Day&Night to Your Blom-Singer® Adhesive Housings

Connect the Day&Night HME to your Blom-Singer Adhesive Housing. To insert the device, hold its outer rim and insert it partially into the lower rim of the previously applied tracheostoma housing. Gently push the device into the applied attachment until it snaps into place. A properly inserted HME cartridge will sit flush and fit snugly inside the housing.

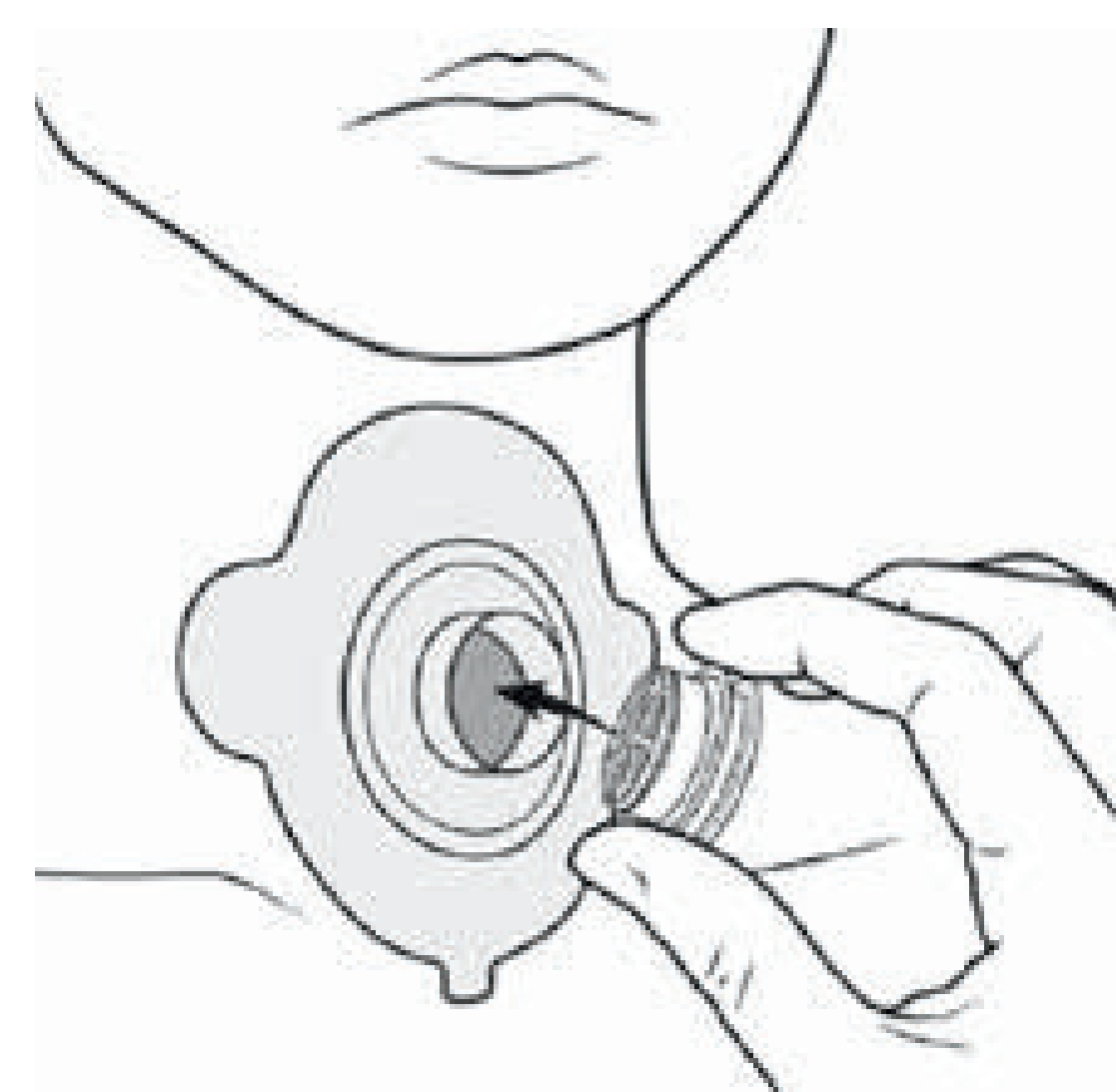

STEP 2:

Test Digital Occlusion


Push with your finger or thumb against the face of the device to talk. Enjoy your HME.



STEP 3:

And at Night...

You can “lock” at night for peace of mind.

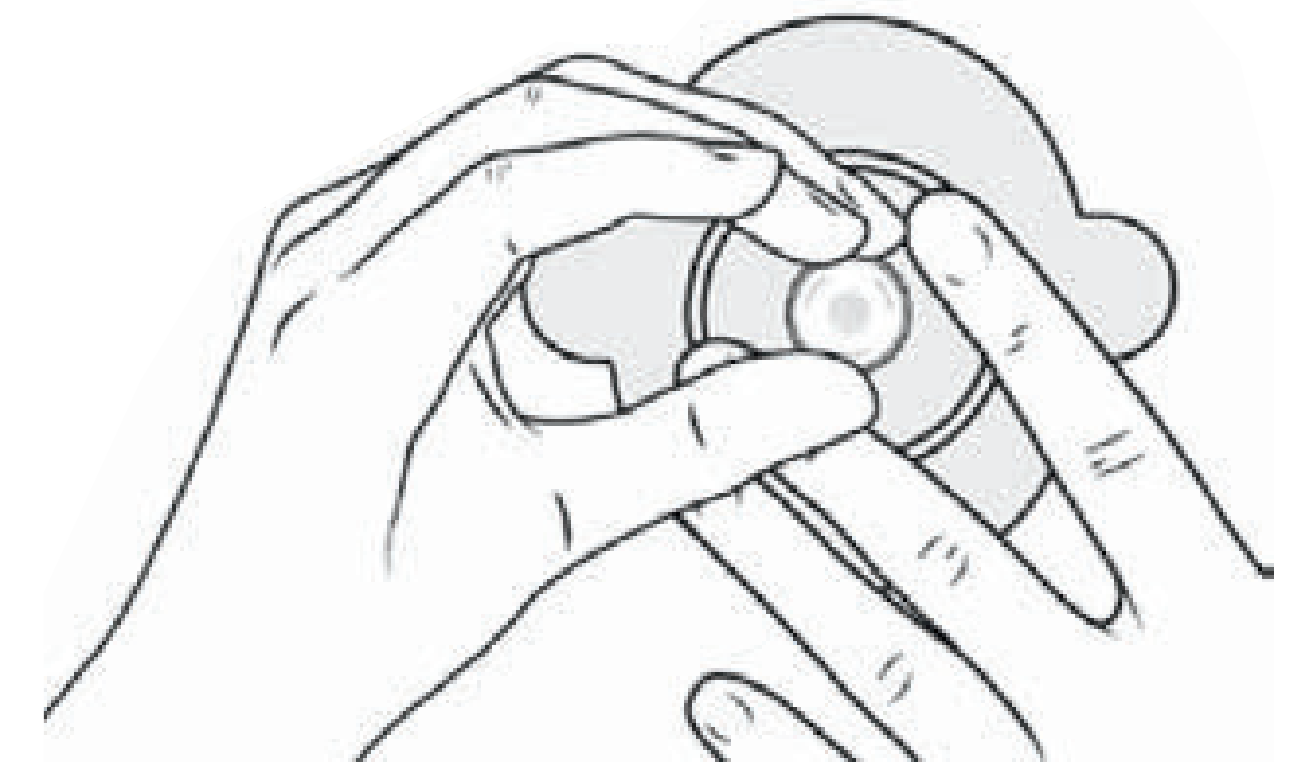

To lock: hold the tracheostoma housing or tube in place with one hand while gently turning the upper lip of the device clockwise with the other hand.

To unlock: hold the tracheostoma housing or tube in place with one hand while gently turning the upper lip of the device counterclockwise with the other hand.

